


Case Report

Prophylaxis for PEP

University Hospital Virgen de Arrixaca
Murcia, Spain

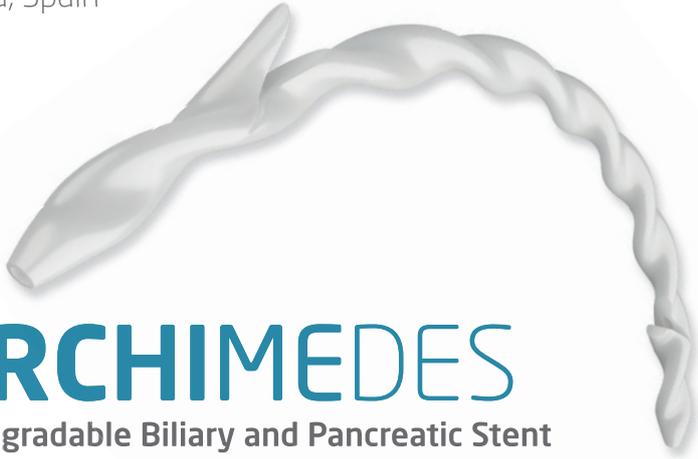

ARCHIMEDES

Biodegradable Biliary and Pancreatic Stent

57 y/o female patient from referred from different hospital facility, with previous repeated episodes of pancreatitis, possibly attributed to pancreas divisum. ERCP performed to evaluate and assess for possible sphincterotomy and pancreatic stent placement.



Image 1. Initial of transpapillary positioning of ARCHIMEDES

Under sedation, endoscopist attempted to assess the biliary tree for diagnostic; however, Papilla of Vater was extremely difficult to cannulate due to the presence of some fibrotic tissue. After 20 minutes, the device was exchanged to a thinner 0.021" sphincterotome system; guide-wire was positioned in the biliary duct and contrast injected showing no alterations. No pancreatic duct (PD) suggestive of pancreas divisum, therefore, also no presence of minor papilla.



Image 2. Implantation of ARCHIMEDES - view with pusher



Image 3. Implantation of ARCHIMEDES - view with pusher



Image 4. Retrieval of pusher, guidewire still in position

Case Report

Prophylaxis for PEP

University Hospital Virgen de Arrixaca
Murcia, Spain

Due to the intense manipulation of the papilla, prophylactic pancreatic stent placement was considered to prevent post-ERCP pancreatitis (PEP). Positioning of the guidewire in the pancreatic duct proved very difficult and an 0.025" guidewire was successfully introduced. An **ARCHIMEDES** biodegradable stent (6 F x 40 mm, fast degradation) was smoothly inserted in the PD in a transpapillary position.



Image 5. Final position of **ARCHIMEDES** biodegradable pancreatic stent

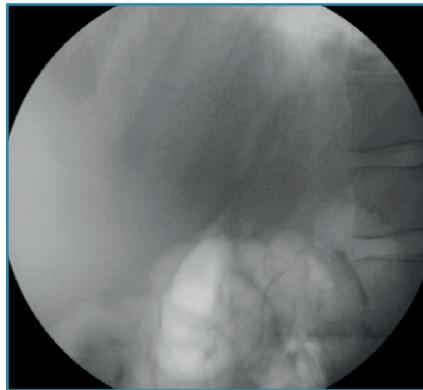


Image 6. Final fluoroscopic view of **ARCHIMEDES** biodegradable stent

The **ARCHIMEDES** biodegradable stent was prepared per the IFU and implanted in the PD, over a 0.025" guidewire, using a 7 F pushing catheter; procedurally similar to the placement of standard plastic stents.

Having the **ARCHIMEDES** biodegradable stent implanted assured that there would be no need for an additional ERCP for a device removal as is the standard with plastic stents.



CE 1434

amginternational
A Q3 Medical Company

amg International GmbH | Boschstraße 16 | D-21423 Winsen | Germany

Phone +49 4171 6905 57-0
Fax +49 4171 6905 57-11

Email info@amggastro.com
Web www.amggastro.com



CREATING VALUE BY HELPING PEOPLE

CCR_BD_GI_BPS-E_003 - Arrixaca_ES - Prophylaxis for PEP - 2023-03
Made in Germany